


Hip Orthosis Traction Device

髖關節矯形牽引器

Model : TD-JX

Features:

- High strength aluminum alloy brace, strong fixed knee joint in extension position
- Three piece design, suitable for people of different statures
- Applicable to:
 1. Fixation of tibial plateau fracture, distal femur fracture, proximal fibula fracture, or ligament injury of knee joint before and after operation
 2. It can replace plaster or splint (it should be used under the guidance of a specialist)
- 高強度鋁合金支條, 強力固定膝關節在伸直位
- 三片式設計, 適用於不同身材人群
- 適用於:
 1. 脛骨平臺骨折、股骨遠端骨折、腓骨近端骨折, 或膝關節各個韌帶損傷手術前後的固定
 2. 可以代替石膏或夾板 (應在專科醫生指導下使用)

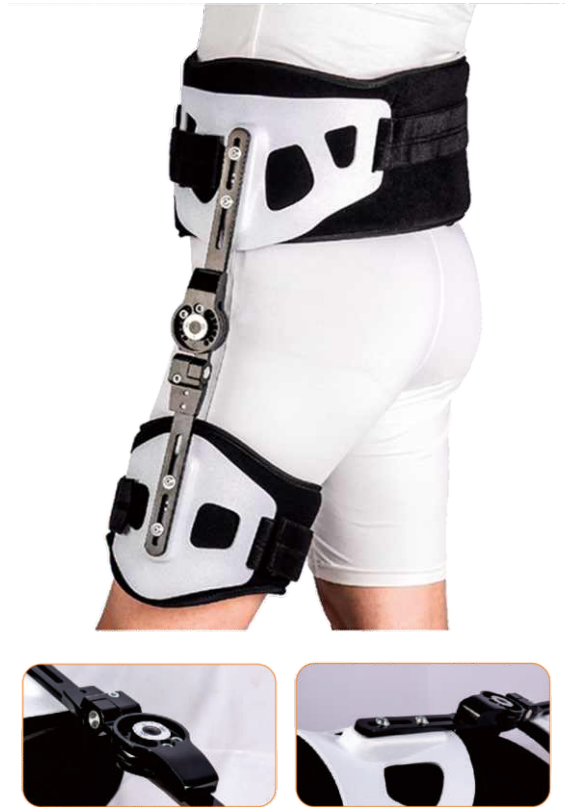

Specifications

Product Model 產品型號	TD-JX
Product Function 產品功能	fixation of knee joint 膝關節固定
Product Material 產品材質	Polyurethane + ABS 聚氨酯+ABS
Fitting 貼合物	100% polyurethane 100%聚氨酯
Specification 規格	free size 均碼
Suitable For 適用於	waist circumference: 80-120cm leg circumference: 50-70cm 適用腰圍: 80-120cm 適用腿圍: 50-70cm
Product Size 產品尺寸	460*270*80mm
Product Weight 產品重量	830g
Package Weight 包裝重量	430*21*270mm
Packing Size 包裝尺寸	830g

Instructions

- ◆ 將固定器上部固定於腰部, 側板位於患肢外側固定器下縛帶固定於大腿。
- ◆ 調整螺絲位置, 可調節固定器長短以適合固定患肢為宜
- ◆ 用配真的螺絲刀(六稜)調節卡盤, 跟據醫囑限定髖關節活動範圍

注意事項:

- ◆ 未成年人、無自理能力者、淺感覺減退者需在監護人幫助下使用
- ◆ 諮詢專業人士正確使用產品
- ◆ 凝血機制障礙、血管脆性增高及精神心理有障礙者不建議使用
- ◆ 此產品佩戴應鬆緊適度, 佩戴後如病情無明顯改善, 請及時諮詢醫生或相關專業人員